

PACS Client Upgrade

Quick Reference Guide

(For Philips Documentation)

DOCUMENT REVISION HISTORY

<i>Revision #:</i>	<i>Revision Date:</i>	<i>Revised by:</i>	<i>Description:</i>
001	Oct 22, 2020	Michael Schultz	Initial Version
002	Oct 23, 2020	Paul Sloman	Updated Version
003	Oct 24, 2020	Lynn Crook	Updated Version

RELATED DOCUMENTS

<i>Title:</i>	<i>Description:</i>
Installing and Upgrading the IntelliSpace PACS 4.4 Client	Philips instructions for Installation and Upgrade of PACS

Contents

Document Classification	2
Purpose of Document.....	2
Uninstalling PACS 4.4.541.5 Instructions	3
Installing IntelliSpace PACS Enterprise or Radiology Clients	4
Launching the PACS Application.....	7
Accessing IntelliSpace PACS Enterprise Client via Internet Explorer.....	9

Document Classification

Internal: information that is generally available to health stakeholders and approved non-health stakeholders

Purpose of Document

This document is intended to be a quick reference guide to help Information Technology resources with the installation/upgrade of Radiology and Enterprise clients for Intellispace PACS 4.4.553.35. This document should be used in conjunction with the above referenced Philips documentation.

Note: A complete uninstallation of the previously installed version of PACS client is a requirement prior to proceeding with the installation of the new version.

Please note this only applies if you installed the Radiology **OR** Enterprise client of PACS on your computer. If you access PACS client via the website, the new upgrade will automatically be installed when you go to the website: <https://pacs.ehealthsask.ca/> by the philips package manager.

This document is a Quick Reference guide **ONLY**. This document is intended for use in conjunction with the Philips Documentation **Installing and Upgrading the IntelliSpace PACS 4.4 Client**. In case of discrepancies, please follow the Philips documentation.

NOTE: If you encounter problems or issues while using this reference sheet or Philips documentation, please call the eHealth Saskatchewan Service Desk at 1-888-316-7446.

Uninstalling PACS 4.4.541.5 Instructions

Please refer to Page 13 of **Installing and Upgrading the IntelliSpace PACS 4.4 Client**


1. Ensure you have logged out of PACS and your browser is closed.
2. Click on the Start Menu.
3. Go to Windows System -> Control Panel -> Programs and Features (Programs>Uninstall a program).
4. Uninstall each of the following. (If installed)
 - a. IntelliSpace PACS Enterprise 4.4.541.xx
 - b. IntelliSpace PACS Radiology 4.4.541.xx
 - c. Philips Medical Systems IntelliSpace PACS Clinical Applications
5. Navigate to C:\Program Files (x86)\Philips\
 - a. Delete the “IntelliSpace PACS *Enterprise*” folder (if it exists)
 - b. Delete the “IntelliSpace PACS *Radiology*” folder (if it exists)
6. Navigate to C:\Users\%username%\AppData\Local\PPM\
 - a. Delete all files and folders within the PPM folder
7. Reboot the computer.

Installing IntelliSpace PACS Enterprise or Radiology Clients


Please refer to Page 25 of **Installing and Upgrading the IntelliSpace PACS 4.4 Client**

Installing the Standalone Client:

1. If you are going through these steps remotely, **please ensure you are first connected to VPN.**
2. Turn off Pop Up Blockers
3. Go to <https://pacs.ehealthsask.ca/installers>
4. Click to the link for “IntelliSpacePACSEnterpriseSetup.exe” or for Radiologist use only “IntelliSpacePACSRadiologySetup.exe”
5. Click “Run” when the window pops up at the bottom:
- Or, you can save the executable file to your local computer and double-click when ready.

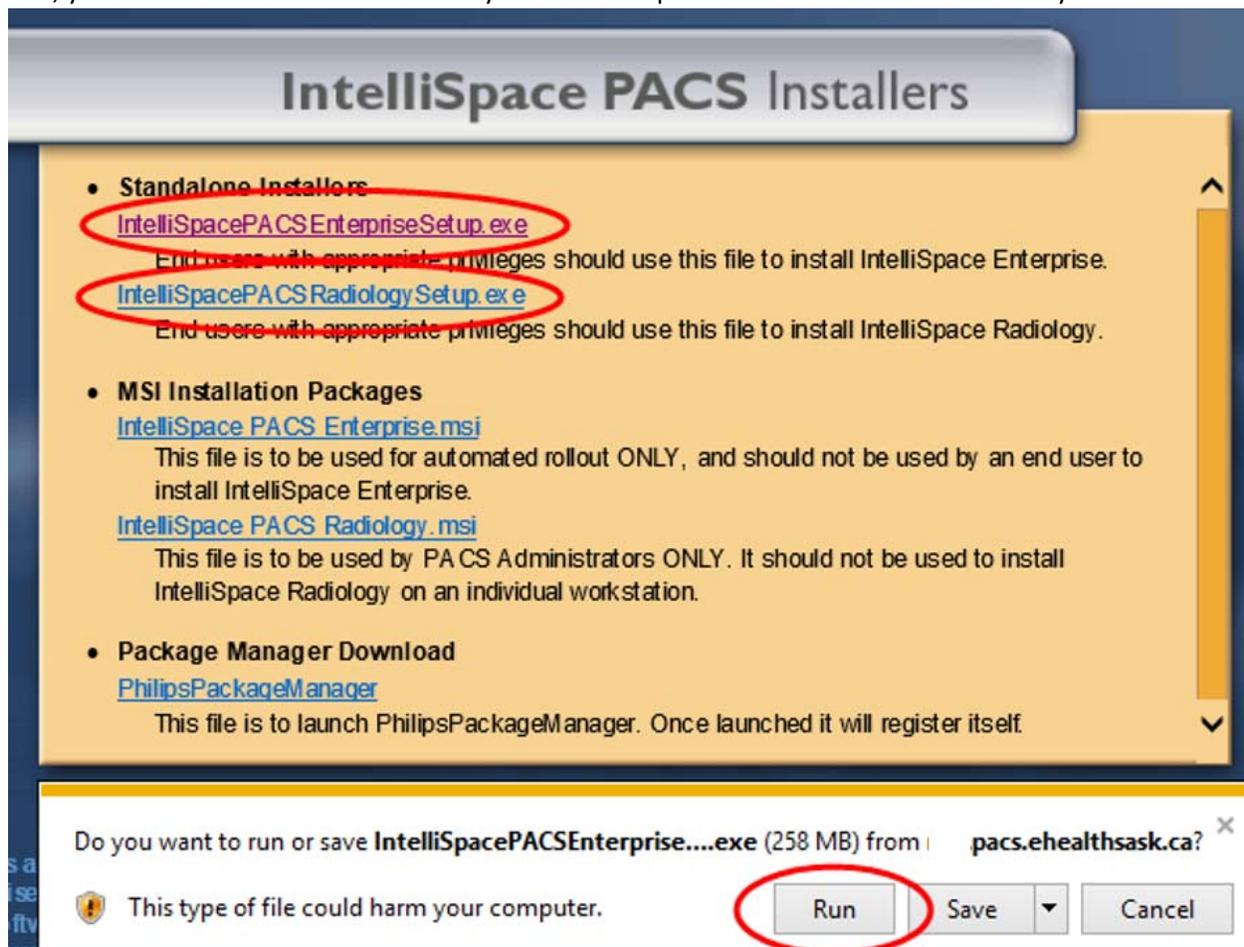


6. Click “Yes” on the “User Account Control” window, if asked

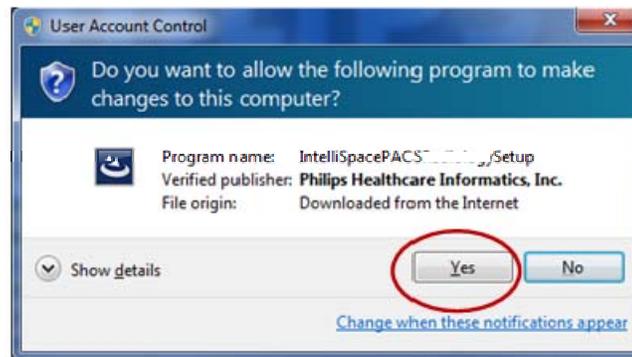


7. When prompted ensure English(United States) is selected, then click “OK”



8. The IntelliSpace PACS Enterprise 4.4.553.35 Install wizard window will pop up. Click “Next”

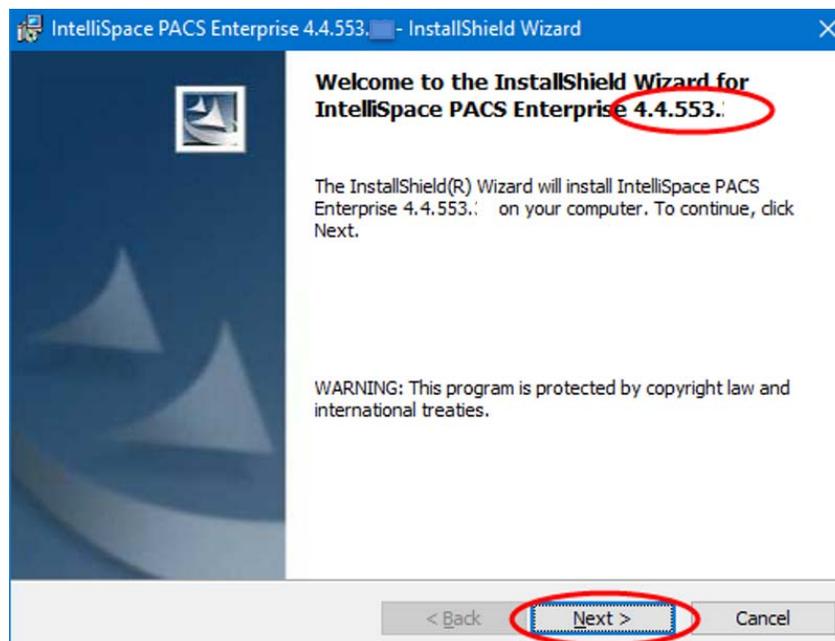


9. Click “Next” on the Custom Setup screen

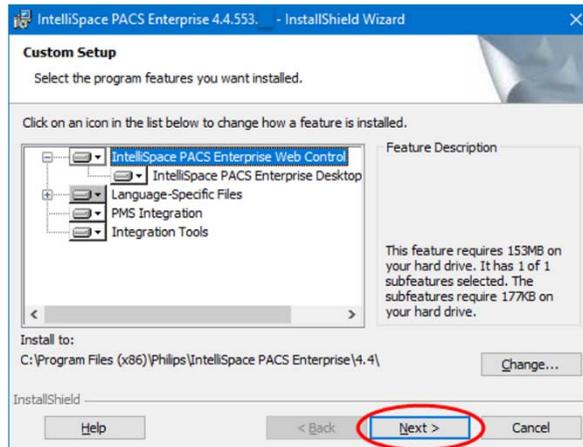


10. Client Configuration screen requires that Server = **pac.s.ehealthsask.ca**.
Then click “Install”


NOTE: If the server prompt does not come up the new version of PACS will **default** to the Server used by the prior version of PACS.

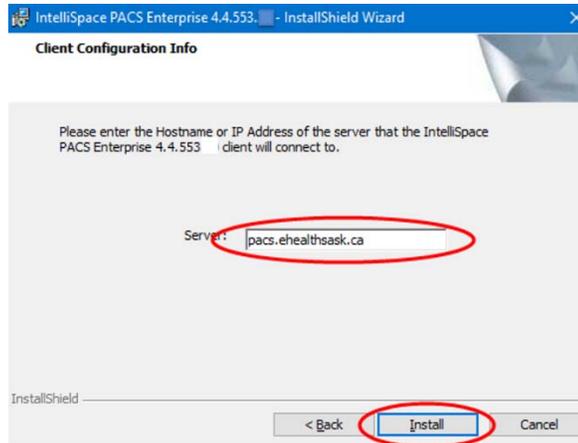

11. A status window will display:

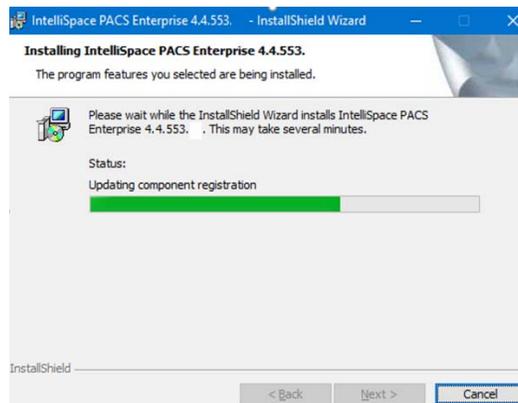

12. When successful, you will get the Install Wizard Completed screen. Click “Finish”

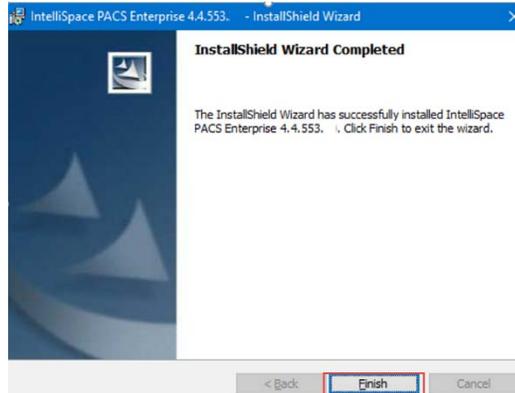


The new version of PACS successfully installed and the application icon will appear on your desktop.

Launching the PACS Application

1. To launch the application, double click the IntelliSpace PACS Enterprise icon on the desktop or start menu:



2. The Philips IntelliSpace Log On Window will appear and display version **4.4.553.35**

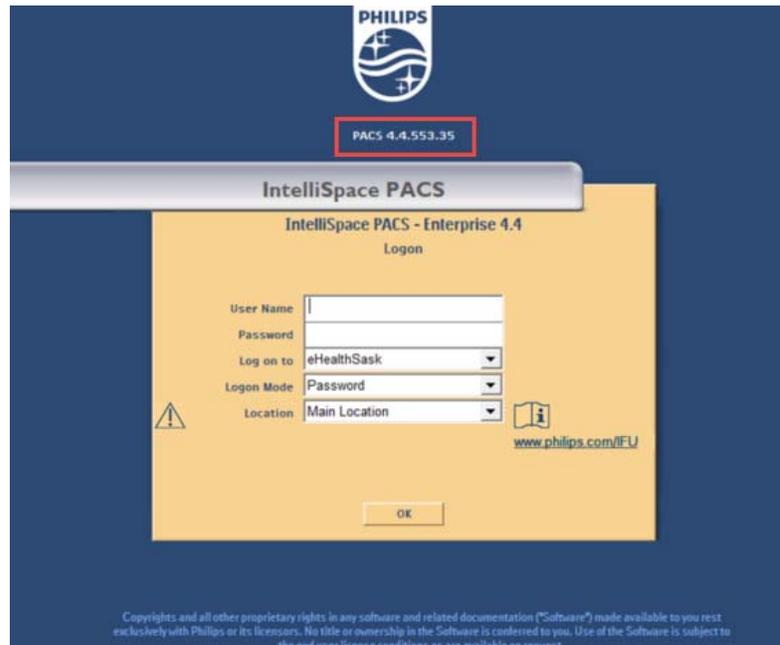


3. Fill out the logon information the same way you have done in the past
- User Name = same as before
 - Password = same as before
 - Log On To = domain where your log on credentials were created
 - Logon Mode = Password
 - Location = Main Location

Accessing IntelliSpace PACS Enterprise Client via Internet Explorer


Please refer to Page 39 of **Installing and Upgrading the IntelliSpace PACS 4.4 Client**

Internet Explorer (IE) 10 or 11 are the supported Web Browsers.


1. Compatibility view should not be enabled
2. From the IE “Tools” navigate to the Compatibility View Settings
3. Internet Explorer > Tools > Compatibility View
4. Ensure that the below options are unchecked

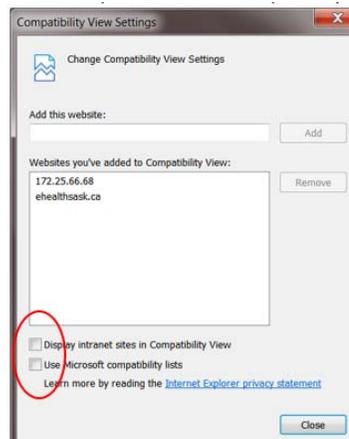


5. To access PACS, ensure you are connected to CNET (Directly or via VPN).
6. Open up Internet Explorer and browse to <https://pacs.ehealthsask.ca>

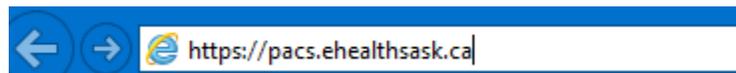

7. On a new installation the Philips Package Manager will launch and display **PACS 4.4.553.35** with a progress bar to show the application is installing. If install does not automatically begin then click on the “**Click here**” link, depending on your browser connection. Please be patient.

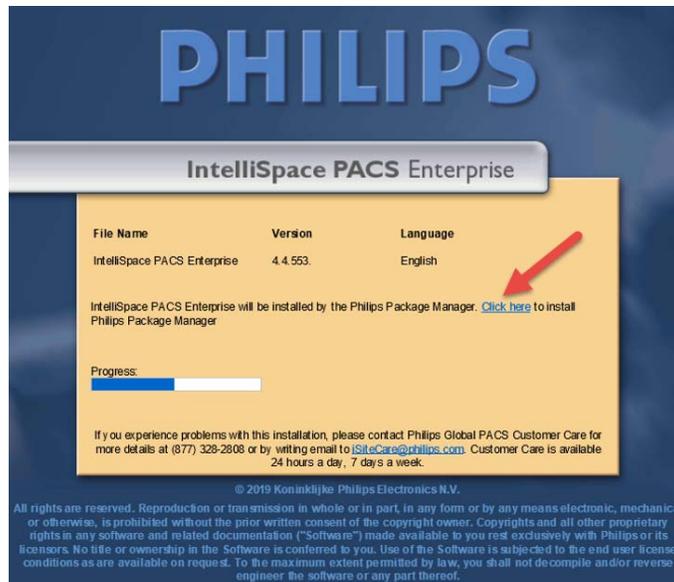

8. Depending on your security permission on the computer, you may have to “allow” the installation:

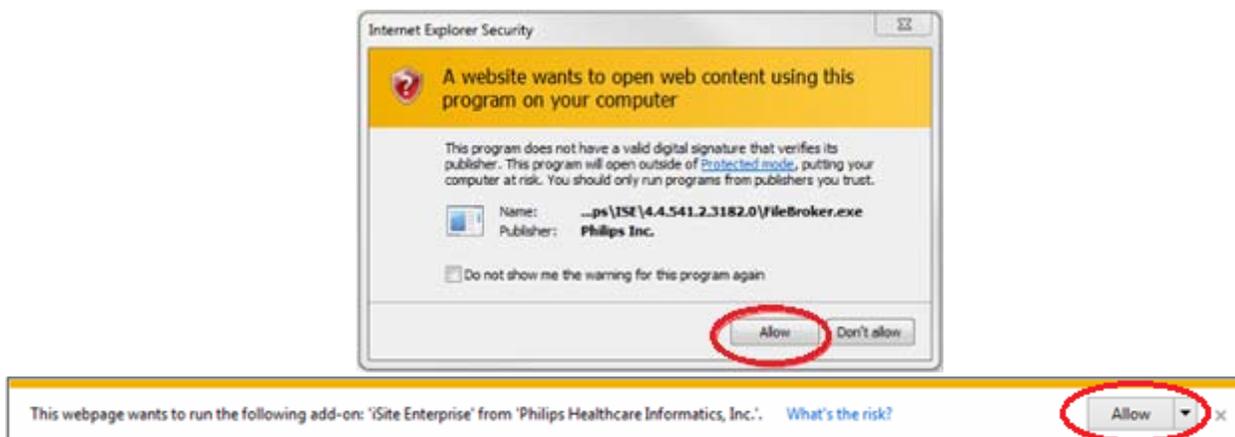

9. After the install is complete. The PACS Logon page will appear as below with PACS version displaying as **PACS 4.4.553.35**

PHILIPS
PACS 4.4.553.35

IntelliSpace PACS

IntelliSpace PACS - Enterprise 4.4
Logon

User Name

Password

Log on to eHealthSask

Logon Mode Password

Location Main Location

www.philips.com/IFU

OK

Copyrights and all other proprietary rights in any software and related documentation ("Software") made available to you rest exclusively with Philips or its licensors. No title or ownership in the Software is conferred to you. Use of the Software is subject to the end user license conditions as are available on request.

10. It is recommended to bookmark this page for easy access.
11. Your log on credentials remain the same
- User Name = same as before
 - Password= same as before
 - Log On To = domain where your log on credentials were created
 - Logon Mode = Password
 - Location = Main